

# seca 676 r

Electronic wheelchair scale with handrail and transport castors and integrated RS232 interface



- Generously sized platform.
- Very stable railing as standing aid.
- Foldable and mobile.
- Space-saving.
- Pre-TARE function with three memory cells.

## seca 676 r:

### Safe weighing – even when standing is difficult.

It is good when one is physically fit and mobile. But unfortunately not everyone is: many people find it increasingly more difficult to walk and stand when they get old. Many suffer from metabolic illnesses such as diabetes and often suffer from overweight as a result. Old and young alike suffer injuries such as, for example, a broken bone; women are strongly impeded in the last months of a pregnancy. Hospitals, surgeries, rehabilitation centres as well as geriatrics and nursing homes must be prepared for all these handicaps. And they are – with the wheelchair scale seca 676 r.



#### Generously sized platform and stable railing.

Thanks to its especially large, non-slip platform and robust rail, the seca 676 r is very versatile. On the one hand, patients can be weighed sitting in a wheelchair or on a chair. On the other hand, people who are fragile or can only walk with great difficulty are able to support themselves.




Thanks to the ramp, a wheelchair rolls easily and comfortably onto the platform.



#### Foldable, mobile and space-saving.

After weighing, the seca 676 r can be folded together in no time at all to save space. The sturdy locking device between rail and platform ensures the scale stands safely also when folded together. The rail also serves as handle when the scale is folded together, meaning that the seca 676 r can be moved and stored away easily and quickly on the smooth-running transport castors.

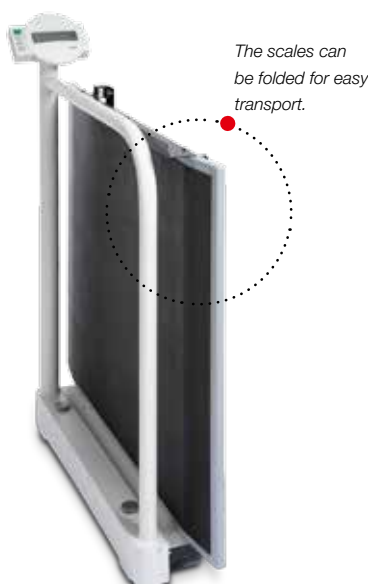

The scales can be folded for easy transport.

#### Intelligent functions for comfortable weighing.

The swivel-mounted display at hip level ensures comfortable operation for doctors and carers. After having been weighed, the patient can first be attended to without any stress as the weight is stored by the HOLD function. The additional weight of all aids used such as chair or wheelchair can be deducted at the touch of a button with the aid of the TARE function. The scale has memory cells for storing three additional weights and the height of three persons. The BMI function provides information on the individual nutritional condition of the weighed person.

## seca 676 r

### Technical Data

- Capacity: 360 kg / 800 lbs
- Graduation: 50 g / 0.1 lbs
- Dimensions (WxHxD): 920 x 1,120 x 1,150 mm / 36.2 x 44.1 x 45.3"
- Dimensions platform (WxHxD): 800 x 55 x 965 mm / 31.5 x 2.2 x 38"
- Weight: 38 kg / 83.8 lbs
- Power supply: power adapter
- Functions: TARE, Pre-TARE, HOLD, Auto-HOLD, BMI, Auto-CLEAR, kg/lbs switch-over, damping, RS232 interface
- Optional: PC interface cable seca 451